


EN
01.2024

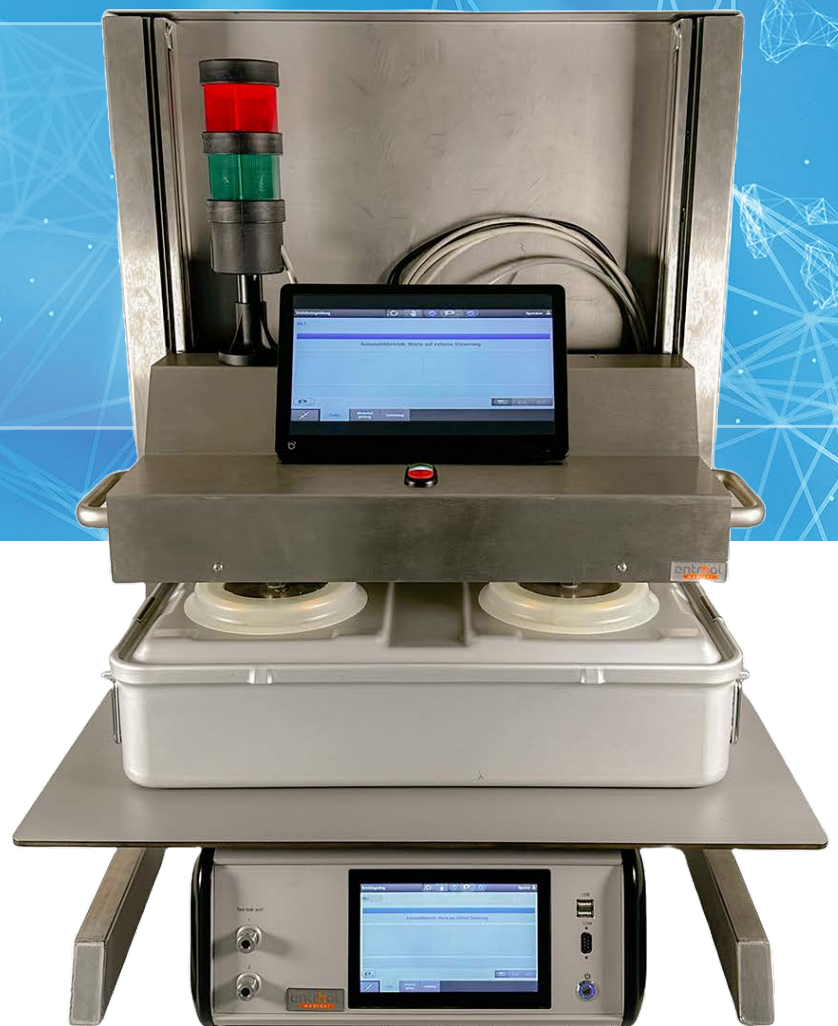

CONTEG
PROFESSIONAL

Container - Testdevice

Security in seconds

www.entrhal-medical.de

entrhal
medical

**Patent
Pending**

- + simple
- + fast
- + reliable
- + reproducible



Fig. Conteg Pro

CONTEG

PROFESSIONAL

The smart testing procedure

Compact container testing device for precise testing of the seal integrity between the Container lid and the base.

The test could be done after packing and before the sterilization process.

Also suitable for quality control in the manufacturing process of rigid sterilization container.

Suitable for most common container brands.

Objective test result in just one minute!!

The strengths of the Conteg container test device

- + Automatic determination of leakage in ml/min
- + Test time only approx. 60 seconds
- + free definition of acceptance levels
- + Storage of measurement data in the device
- + Network capable
- + Compact design
- + Clear display of test results with green (pass) or red (fail) signal lights
- + Extended detailed view on separate 12 inch screen (optionally also available in the base version)
- + Service module slide-in system for easy maintenance and calibration of the measuring unit
- + Individual adaptation of different container brands possible


Available versions

CONTEG PROFESSIONAL

Version: Pro

Mechanical covering of container openings by push and pull device

Article-No: 8881176

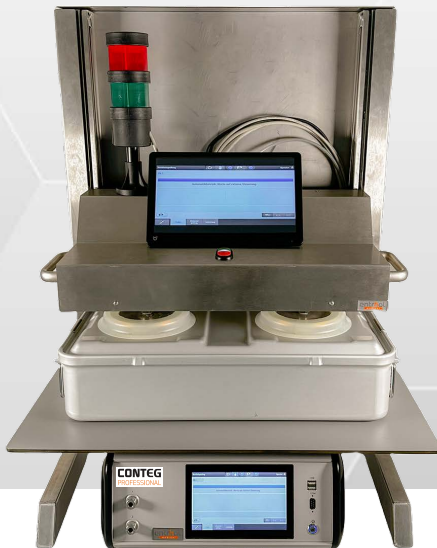

- + Automatic determination of leakage in ml/min
- + Test time only approx. 60 seconds
- + free definition of acceptance levels
- + Storage of measurement data in the device
- + Network capable
- + Compact design
- + Always correctly positioned thanks to standardised guidance
- + Clear display of test results with green (pass) or red (fail) signal lights
- + Extended detailed view on separate 12 inch screen
- + Service module slide-in system for easy maintenance and calibration of the measuring unit
- + Individual adaptation of different container brands possible
- + Testing of different container sizes (S,M,L)

CONTEG PROFESSIONAL

Version: Base

Manual covering of the container openings by attaching the sealing bells

Article-No: 8881175



- + Automatic determination of leakage in ml/min
- + Test time only approx. 60 seconds
- + free definition of passed acceptance levels
- + Storage of measurement data in the device
- + Network capable
- + Compact design
- + Always correctly positioned thanks to standardised guidance
- + Clear display of test results with green (pass) or red (fail) signal lights
- + Optionally available with separate 12 inch display
- + Service module slide-in system for easy maintenance and calibration of the measuring unit
- + Individual adaptation of different containers possible
- + Testing of different container sizes (S,M,L)

CONTEG

PROFESSIONAL

Features

Dimensions	Version Pro	Version Base
WxDxH	69 x 45 x 102 cm	
Housing	Stainless steel	Stainless steel
Mains voltage	230 V	230 V
Mains voltage	50 Hz	50 Hz
Power supply	24 VDC (+10 / -5 %) 5 A	24 VDC (+10 / -5 %) 5 A
Test medium	Compressed air (oil and water free according to ISO 8573-1 Class 3)	Compressed air (oil and water free according to ISO 8573-1 Class 3)

Sealing devices + accessories

Test equipment

Instrument cleaning

Disposable instruments + sets

Entrhal Medical GmbH | Fax: +49 (0) 2834 942 48-46
Phone: +49 (0) 2834 942 48-0 | info@entrhal-medical.de

www.entrhal-medical.de

entrhal
medical