

Symmetry Surgical's Quad-Lock® container system now features **NEW** sterilization modalities with a patented, advanced design engineered to ensure dependable lid attachment.

- Four latches secure the lid
- Validated and FDA cleared for pre-vacuum steam, ethylene oxide and low-temperature hydrogen peroxide sterilization modalities
- Half-size containers and extra-deep full-size containers
- Ergonomic silicone handles with 90° stops, secondary security lids, and multiple color options



Quad-Lock® Container System


Quad-Lock® Quick Reference

Container System by Sterilization Modality



Quad-Lock® Sterilization Container

Baskets

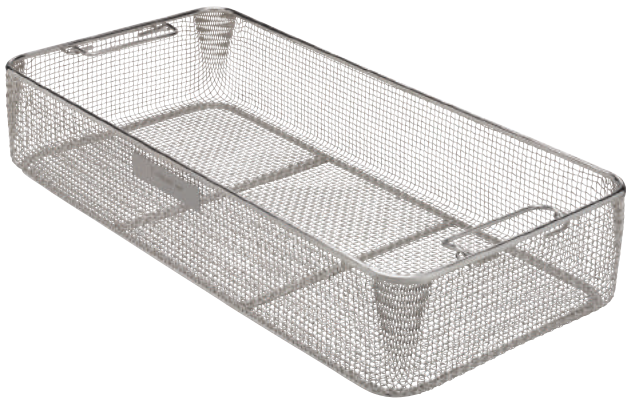

- Sturdy lightweight holders for your instrument sets
- Available in a number of sizes
- Designed with handles that feature a 90° stop for easy access
- Handles can also slide neatly out of the way, for easy instrument access

Please consult product labels and inserts for indications, contraindications, hazards, warnings, cautions and instructions for use.

	Pre-Vacuum/ Dynamic Air	Ethylene Oxide	Low-Temperature Hydrogen Peroxide
Basket	X	X	X
Filtered Bottom	X	X	X
Solid Bottom	X		
Wire Base	X	X	X
Grey Container Lid	X	X	X
Colored Container Lid	X		
Protective Mats	X		
Security Lids	X	X	X

Visit [symmetrysurgical.com](https://www.symmetrysurgical.com) or contact your local Sales Representative for more information.

Baskets

CODE	DESCRIPTION	HEIGHT
50-8500	Full Length	30mm (1 inch)
50-8501	Full Length	50mm (2 inch)
50-8502	Full Length	70mm (3 inch)
50-8503	Full Length	100mm (4 inch)
50-8504	Half Length	30mm (1 inch)
50-8505	Half Length	50mm (2 inch)
50-8506	Half Length	70mm (3 inch)
50-8507	Half Length	100mm (4 inch)
50-8508	Three-Quarter Length	30mm (1 inch)
50-8509	Three-Quarter Length	50mm (2 inch)
50-8510	Three-Quarter Length	70mm (3 inch)
50-8511	Three-Quarter Length	100mm (4 inch)

Quad-Lock® Sterilization Container Base



- Number of heights available to accommodate even your largest instrument set
- The oversized ergonomic handles have a 90° stop to protect fingers during transportation
- The container base has specially-designed feet for secure stacking and standard slots for large or small indicator cards or aluminum name plates
- All Filter Container Bases come fully equipped with Filter Retention Discs

Selecting the Container Height

- Add the height of the selected basket, plus the height of the wire base, plus an additional 50mm (2 Inches) of required space
- For example: 2 Inch Basket + 1 Inch Wire Base + 2 Inches = 5 Inch Container Base
- See product Instructions For Use for further details

Container Base Dimensions

BASE	LENGTH	WIDTH
Full Length	560mm (22 inches)	275mm (11 inches)
Three-Quarter	440mm (17.5 inches)	270mm (10.5 inches)
Half-Size	270mm (10.5 inches)	270mm (10.5 inches)

Visit [symmetrysurgical.com](https://www.symmetrysurgical.com) or contact your local Sales Representative for more information.

Filter Bottom Container Base

Filter Bottom Container Bases must be used for the ethylene oxide and low temperature hydrogen peroxide gas sterilization modalities.

CODE	DESCRIPTION	HEIGHT
50-8726	Full Length	100mm (4 inch)
50-8727	Full Length	135mm (5 inch)
50-8728	Full Length	150mm (6 inch)
50-8729	Full Length	200mm (8 inch)
50-8730	Full Length	260mm (10 inch)
50-8731	Half Length	100mm (4 inch)
50-8732	Half Length	135mm (5 inch)
50-8733	Half Length	150mm (6 inch)
50-8734	Half Length	200mm (8 inch)
50-8735	Half Length	260mm (10 inch)
50-8736	Three-Quarter Length	100mm (4 inch)
50-8737	Three-Quarter Length	135mm (5 inch)
50-8738	Three-Quarter Length	150mm (6 inch)

Solid Bottom Container Base

CODE	DESCRIPTION	HEIGHT
50-8739	Full Length	100mm (4 inch)
50-8740	Full Length	135mm (5 inch)
50-8741	Full Length	150mm (6 inch)
50-8742	Full Length	200mm (8 inch)
50-8743	Full Length	260mm (10 inch)
50-8744	Half Length	100mm (4 inch)
50-8745	Half Length	135mm (5 inch)
50-8746	Half Length	150mm (6 inch)
50-8747	Half Length	200mm (8 inch)
50-8748	Half Length	260mm (10 inch)
50-8749	Three-Quarter Length	100mm (4 inch)
50-8750	Three-Quarter Length	135mm (5 inch)
50-8751	Three-Quarter Length	150mm (6 inch)

Please consult product labels and inserts for indications, contraindications, hazards, warnings, cautions and instructions for use.

Visit [symmetrysurgical.com](https://www.symmetrysurgical.com) or contact your local Sales Representative for more information.

Quad-Lock® Sterilization Container


Wire Base



- Designed to prevent instruments, protective mats and baskets from touching the bottom of the container, for proper sterilization and air circulation purposes
- A Wire Base must be used with every Container Base that has filters on the bottom

Quad-Lock® Sterilization Container

Protective Mats

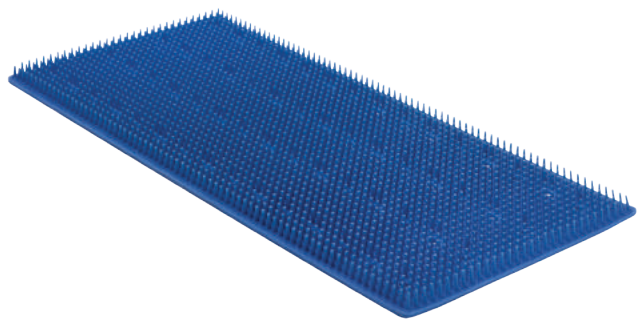

- Helps protect the most delicate instruments
- Pair a Protective Mat with a Wire Base, and you can place large objects directly in the Container Base without using a Basket
- Protective mats are not compatible with all modalities of low temperature hydrogen peroxide sterilization. See Product Instructions For Use for further details.

Please consult product labels and inserts for indications, contraindications, hazards, warnings, cautions and instructions for use.

Wire Bases

CODE	LENGTH	HEIGHT
50-8512	Half Length	26mm (1 Inch)
50-8513	Three-Quarter Length	26mm (1 Inch)
50-8514	Full Length	26mm (1 Inch)

Protective Mats

CODE	LENGTH
50-8550	Half Length
50-8551	Three-Quarter Length
50-8552	Full Length


Quad-Lock® Sterilization Container Filter Lids



- Advanced, patented Quad-Lock® design for reliable lid attachment
- Four latch lids available in six colors, feature a specially molded design to avoid warping without adding weight
- Deep-seated grooves for secure stacking

Lid colors available include: Red, Yellow, Green, Blue, Grey and Black

(Colored lids may be used for steam sterilization ONLY. Grey filter lids compatible with all sterilization modalities)

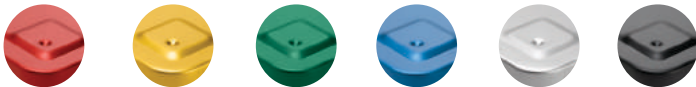

Quad-Lock® Sterilization Container Security Lids



- Provide an additional layer of protection for the filter
- Self attaching 'drop on' latches, deep seated grooves for secure stacking
- Made of Anodized Aluminum for lightweight durability

Lid colors available include: Red, Yellow, Green, Blue, Grey and Black



Please consult product labels and inserts for indications, contraindications, hazards, warnings, cautions and instructions for use.

Filter Lids

Colored filter lids may be used with steam sterilization ONLY. Grey filter lids compatible with all sterilization modalities.

CODE	DESCRIPTION	COLOR
50-8944	Full Length	Grey
50-8945	Full Length	Yellow
50-8946	Full Length	Green
50-8947	Full Length	Blue
50-8948	Full Length	Red
50-8949	Full Length	Black
50-8950	Half Length	Grey
50-8951	Half Length	Yellow
50-8952	Half Length	Green
50-8953	Half Length	Blue
50-8954	Half Length	Red
50-8955	Half Length	Black
50-8956	Three-Quarter Length	Grey
50-8957	Three-Quarter Length	Yellow
50-8958	Three-Quarter Length	Green
50-8959	Three-Quarter Length	Blue
50-8960	Three-Quarter Length	Red
50-8961	Three-Quarter Length	Black

Security Lids

CODE	DESCRIPTION	COLOR
50-8869	Full Length	Grey
50-8870	Full Length	Yellow
50-8871	Full Length	Green
50-8872	Full Length	Blue
50-8873	Full Length	Red
50-8874	Full Length	Black
50-8875	Half Length	Grey
50-8876	Half Length	Yellow
50-8877	Half Length	Green
50-8878	Half Length	Blue
50-8879	Half Length	Red
50-8880	Half Length	Black
50-8881	Three-Quarter Length	Grey
50-8882	Three-Quarter Length	Yellow
50-8883	Three-Quarter Length	Green
50-8884	Three-Quarter Length	Blue
50-8885	Three-Quarter Length	Red
50-8886	Three-Quarter Length	Black

Quad-Lock® Sterilization Container

Filters, Security Seals, Labels and Ergonomic Hand Grips



Ergonomic Hand Grip

- Great comfort grip
- Available in six colors (shown below)
- Provide protection from the risks associated with high temperatures following sterilization

Aluminum Name Plates

- Etched to your specific needs, up to 45 characters
- Made of Anodized Aluminum for durability with lightweight construction

Colors available include: Red, Yellow, Green, Blue, Grey and Black



Single-Use Products

- Clearly-labeled round filters, small and large indicator cards all include a sterility indicator

Visit symmetrysurgical.com or contact your local Sales Representative for more information.

Ergonomic Silicone Hand Grips (Pair)

CODE	DESCRIPTION	COLOR
50-9012	Ergonomic Silicone Hand Grips	Red
50-9013	Ergonomic Silicone Hand Grips	Blue
50-9014	Ergonomic Silicone Hand Grips	Green
50-9015	Ergonomic Silicone Hand Grips	Yellow
50-9016	Ergonomic Silicone Hand Grips	Grey
50-9017	Ergonomic Silicone Hand Grips	Black

Aluminum Name Plates (1 each)

CODE	DESCRIPTION	COLOR
50-9000	Aluminum Name Plate	Red
50-9001	Aluminum Name Plate	Blue
50-9002	Aluminum Name Plate	Green
50-9003	Aluminum Name Plate	Yellow
50-9004	Aluminum Name Plate	Grey
50-9005	Aluminum Name Plate	Black


Single-Use Products

CODE	DESCRIPTION	QTY
50-9006	Paper Label with Indicator Small	250
50-9007	Paper Label with Indicator Large	250
50-9008	Security Seal Single-Use Only – Blue	100
50-9010	Filter – Single-Use Only	1000
50-9044	Filter Retention Disc	1

Please consult product labels and inserts for indications, contraindications, hazards, warnings, cautions and instructions for use.

Quad-Lock® Sterilization Container

Accessories

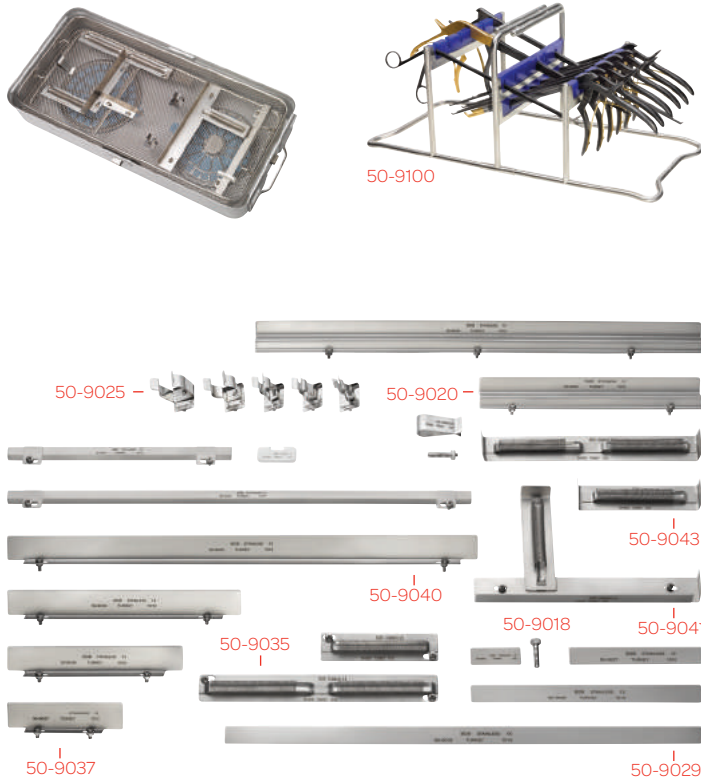

- Customize your containers to fit your exact needs and preferences
- Customize the arrangement of your tray contents for convenience and the protection of your medical equipment
- Accessories are easily installed with an included nut and bolt

Please consult product labels and inserts for indications, contraindications, hazards, warnings, cautions and instructions for use.

Visit [symmetrysurgical.com](https://www.symmetrysurgical.com) or contact your local Sales Representative for more information.

Accessories

CODE	DESCRIPTION
50-9018	Holding Pin for Divider – 1 each (Use with 509026 to 509029)
50-9019	Holding Pin for Instruments 26mm (1 inch) (1 each)
50-9020	Basket Handle Tag (1 each)
50-9021	Spring Clip 7mm–12mm ($\frac{1}{4}$ inch– $\frac{1}{2}$ inch)
50-9022	Spring Clip 12mm–16mm ($\frac{1}{2}$ inch– $\frac{3}{4}$ inch)
50-9023	Spring Clip 16mm–26mm ($\frac{3}{4}$ inch–1 inch)
50-9024	Spring Clip 26mm–36mm (1 inch–1 $\frac{1}{2}$ inch)
50-9025	Spring Clip 36mm–45mm (1 $\frac{1}{2}$ inch–1 $\frac{3}{4}$ inch)
50-9026	Divider 50mm x 20mm (2 inch x $\frac{3}{4}$ inch) (Use with two 509018)
50-9027	Divider 130mm x 20mm (5 inch x $\frac{3}{4}$ inch) (Use with two 509018)
50-9028	Divider 225mm x 20mm (9 inch x $\frac{3}{4}$ inch) (Use with two 509018)
50-9029	Divider 460mm x 20mm (18 inch x $\frac{3}{4}$ inch) (Use with two 509018)
50-9030	Polymer Spacer – 1 each (Use with 509031 and 509032)
50-9031	Sliding Rail for Polymer Spacer 230mm (9 inch) (Use with 509030)
50-9032	Sliding Rail for Polymer Spacer 470mm (18 $\frac{1}{2}$ inch) (Use with 509030)
50-9033	Limiting Bar 230mm (9 inch) (Screws to Basket)
50-9034	Limiting Bar 470mm (18 $\frac{1}{2}$ inch) (Screws to Basket)
50-9035	Spiral Holding Device - Double 230mm (9 inch) (Screws to Basket)
50-9036	Spiral Holding Device - Single 130mm (5 inch) (Screws to Basket)
50-9037	Limiting Bar Right Angle 110mm (4 $\frac{1}{2}$ inch) (Screws to Basket)
50-9038	Limiting Bar Right Angle 170mm (6 $\frac{1}{2}$ inch) (Screws to Basket)
50-9039	Limiting Bar Right Angle 230mm (9 inch) (Screws to Basket)
50-9040	Limiting Bar Right Angle 470mm (18 $\frac{1}{2}$ inch) (Screws to Basket)
50-9041	Spiral Divider with Limiting Bar 130mm x 230mm (5 inch x 9 inch)
50-9042	Spiral Divider - Double 230mm (9 inch) (Use with 509018)
50-9043	Spiral Divider - Single 130mm (5 inch) (Use with 509018)
50-9100	Universal Rongeur Instrument Holder (Full Size Basket)
50-9101	Universal Rongeur Instrument Holder (Half Size Basket)

How to choose the right container:

1. Determine the appropriate size basket to use based on the number, size and weight of the instruments to be sterilized.

NOTE: No instruments should extend above the inner edge of the basket.

2. Based on the sterilization systems used at your facility, determine if a **Solid Bottom Container Base** or **Filter Bottom Container Base** is appropriate.

NOTE: When selecting the container height, add the height of the selected Basket, plus the height of the Wire Base (if used), plus an additional 50mm (2 Inches) of required space.

For example: 2 Inch Basket + 1 Inch Wire Base + 2 Inches = 5 Inch Container Base.

See Product Instructions For Use for further details.

3. Select a **Wire Base** the same length as the basket when using a Filter Base.
4. Based on the sterilization systems used at your facility, choose a **Protective Mat** in the same length as your basket or upgrade to Accessory Instrument Organizers appropriate for your specific instrument set. See Product Instructions For Use for further details.
5. Add a **Filter Lid** in the color you prefer for each instrument set. Filter Lids are available in red, yellow, green, blue, grey and black. Colored lids may be used with steam sterilization ONLY. Grey lids can be used with all sterilization modalities.
6. Select a **Security Lid** in the color you prefer for each instrument set. Security Lids are available in red, yellow, green, blue, grey and black.
7. Consider **Ergonomic Hand Grips** in the color you prefer for each instrument set. Ergonomic Hand Grips are available in red, yellow, green, blue, grey and black.
8. Select **Aluminum Name Plates** in the color you prefer for each instrument set. Aluminum Name Plates are available in red, yellow, green, blue, grey and black. Make note of what should be etched on each Aluminum Name Plate, up to a 45 character limit.

9. DON'T FORGET:

- Quad-Lock® Single-Use Filters
- Quad-Lock® Single-Use Paper Label with Indicator
- Quad-Lock® Single-Use Security Seals



3034 Owen Drive | Antioch, TN 37013 | 615.964.5558 | 800.251.3000
customerservice@symmetrysurgical.com | [symmetrysurgical.com](https://www.symmetrysurgical.com)

©2016 Symmetry Surgical, Inc. All rights reserved. Cat No. SM 0019 Rev 1